

## CORIO CD-B33 Heating circulator with stainless steel open bath


The open heating bath circulators of the new CORIO series feature high-quality bath tanks made of stainless steel.



### Product features

- Models for internal applications
- Bright, white, easy to read display
- Very quiet
- Easy switching between internal and external circulation
- High-quality bath tanks made of stainless steel
- Integrated drain screw (except B5)
- Integrated drain screw

### Heat-up time

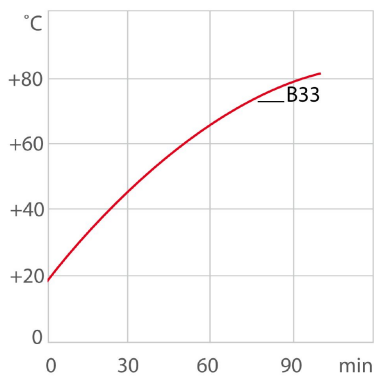

Medium: Water

### Performance values

230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F)	
Heating capacity kW	2
Viscosity max. cSt	50
Pump capacity flow rate l/min	15
Pump capacity flow pressure bar	0.35
Power A	10

### Technical data

Available voltage versions		Bath	
Order No.	9 012 433	Bath tank	Stainless steel
Available voltage versions:		Bath cover	accessory
9012433.01	100V/50-60Hz (Nema N5-15 Plug)	Usable bath opening cm (W x L / D)	66 x 32 / 15
9012433.02	115V/60Hz (Nema N5-15 Plug)		
9012433.03	230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F)		

9012433.03.chn	230V/50Hz (CN Plug)
9012433.04	230V/50Hz (UK Plug Type BS1363A)
9012433.05	230V/50Hz (CH Plug Type SEV 1011)
9012433.13	230V/60Hz (Schuko Plug - CEE 7/4 Plug Type F)

Cooling	
Cooling coil	accessory

Other	
Classification	Classification III (FL)
IP Code	IP 21
Pump function	Pressure Pump
Pump type	Immersion Pump

Electronics	
Interfaces	USB
Temperature control	PID1
Absolute temperature calibration	1 Point Calibration
Temperature display	LED
Temperature setting	Keypad

Dimensions and volumes	
Weight kg	21
Barbed fittings inner diameter	8/12 mm
Dimensions cm (W x L x H)	91 x 36 x 43
Filling volume l	26 ... 39
Pump connections	M16x1 male

Temperature values	
Setting the resolution of the temperature display °C	0.1
Working temperature range °C	+20 ... +150
Temperature stability °C	±0.03
Ambient temperature °C	+5 ... +40
Temperature display resolution °C	0.01 ... 0.1

**Included in delivery**  
 2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections M16x1 male)

For applications near or below ambient temperature: use a cooling coil or JULABO immersion cooler.

### All Benefits



**100% Checked.**  
 100% testing. 100% quality. Each JULABO Circulator undergoes thorough quality testing before leaving the factory.



**ATC.**  
 Absolute Temperature Calibration, 1-point calibration (CD).




**Green technology.**  
 Development consistently applied environmentally friendly materials and technologies.



**Internal. External.**  
 The pump is controlled via a lever located directly below the display. Easily change between internal and external circulation.



**JULABO. Quality.**  
 Highest standards of quality for a long product life.

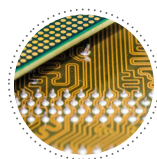

**Modern. Reliable.**  
 High-grade components in every JULABO Circulator – platinum sensors, proven motor technology, CAN-Bus communication and much more.



**Satisfied customers.**  
11 subsidiaries and more than 100 partners worldwide guarantee fast and qualified JULABO support.



**Services 24/7.**  
Around the clock availability. You can find suitable accessories, data sheets, manuals, case studies, and more at [www.julabo.com](http://www.julabo.com).



**Stable.**  
Rubber feet allow for a secured footing of your CORIO to prevent damage to your laboratory equipment.



**Quick start.**  
Individual JULABO consultation and comprehensive manuals at your disposal.



**Timer. Integrated.**  
CORIO circulators include an integrated timer function. When the set time has elapsed, a signal sounds and the device switches off. Setting range: 0 ... 999 minutes.



**Brilliant.**  
Very bright display makes it easy to read even from a distance.



**Connection. Easy.**  
Inclined pump connections (M16×1) facilitate the connection of applications. Each unit includes 2 barbed fittings of 8/12 mm diameter each.



**Everything at the front.**  
All operating controls and safety functions are accessed easily and comfortably from the front.



**Exact.**  
You can rely on it. PID1 control and 'Active Cooling Control' make the new CORIO precise and perfect.



**Locked in.**  
The lockable power plug guarantees safe connection. More process safety.



**Switch on. And off you go.**  
Intelligent operating concept. Ready for operation with just a few quick and easy steps.



**Early warning system for low liquid level**  
Maximum safety for applications, optical and audible alarm, allows user to refill bath fluid before the unit shuts down



**Connectivity.**  
Remote control made easy. CORIO CD circulators feature a USB connection.